


STRIPSPIN™ 12

Mini Centrifuge

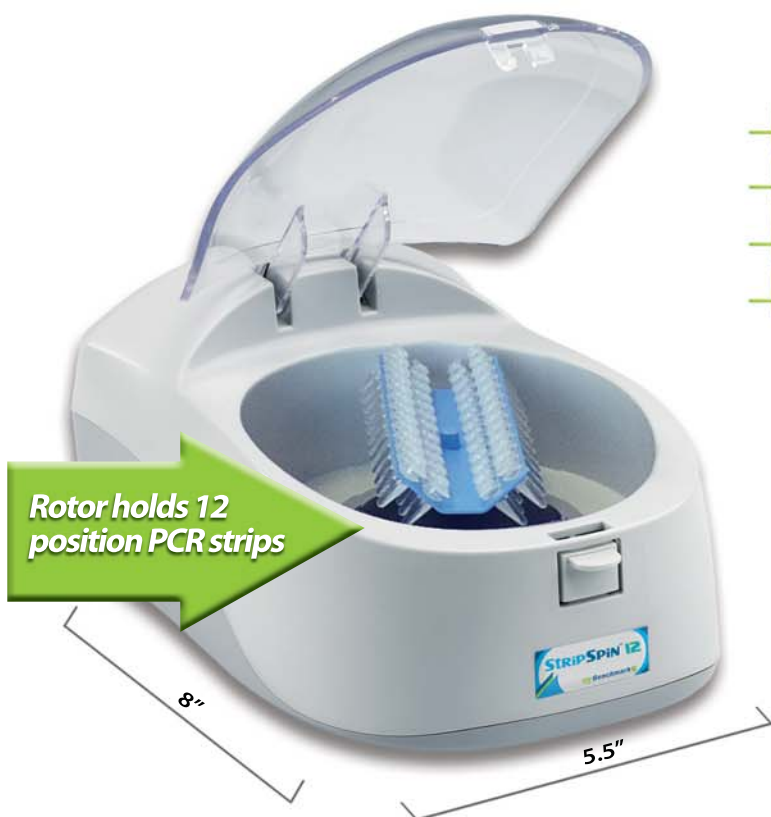


- + Unique rotor holds 12 position PCR strips
- + "Mini" sized, less than 6 inches wide
- + Near silent operation
- + Starts & stops with opening/closing of the lid

Technical Data

Speed:	5,500rpm
Max. g-force:	2,000 xg
Capacity:	4 x 12 position PCR strips 4 x 8 position PCR strips 48 x 0.2ml PCR tubes
Dimensions:	5.5 x 7.9 x 4.4 in 14 x 20 x 11.2 cm
Weight:	2.1 lbs / 5 kg
Electrical:	100 to 240V, 50-60Hz

The new StripSpin 12 Mini Centrifuge from Benchmark is the only mini centrifuge designed for centrifugation of 12 position PCR strips. With a fixed speed that produces 2,000 xg, this centrifuge is perfect for quick spin downs. Simply close the lid and the unit quickly ramps up to 5500 rpm. Open the lid, and the rotor quickly decelerates for removal of samples.

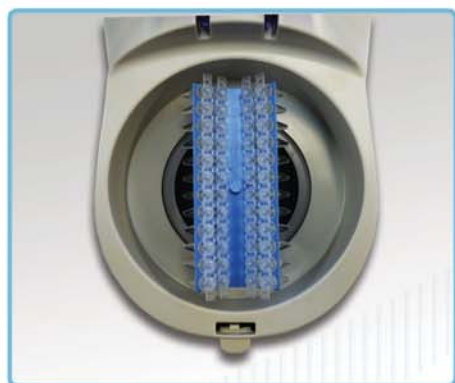

Holds up to 4 x PCR strips or 48 x 0.2ml tubes

Ordering Information

C1248*	StripSpin Mini Centrifuge, with rotor for up to 4 x PCR strips (8 or 12 position) or 48 x 0.2ml tubes
C1008-ROT	Optional rotor, 8 x 1.5/2.0ml
C1008-A5	Optional 0.5ml tube adapters for C1008-ROT, 8/pk.

*To order in 230V, add (-E) to the item number.

Benchmark
Scientific



Benchmark Scientific Inc.
PO Box 709
Edison, NJ 08818
Phone: 908-769-5555
Fax: 908-222-1864

Web: www.BenchmarkScientific.com
Email: Info@BenchmarkScientific.com